



AQA Chemistry GCSE

Required Practical 2

Neutralisation [Foundation Tier]

Methods taken from the AQA Required Practical Handbook



Neutralisation

This method allows the determination of the reacting volumes of solutions of a strong acid and a strong alkali by titration.

Aim

Carry out an investigation to find the volume of dilute sulfuric acid needed to neutralise a known volume of sodium hydroxide solution.

Equipment list

- 25 cm³ volumetric pipette
- Pipette filler
- 50 cm³ burette
- 250 cm³ conical flask
- Small funnel
- Clamp stand and clamp
- White tile
- 0.1 M sodium hydroxide solution
- 0.08 M sulfuric acid
- Phenolphthalein indicator


Method

1. Use the pipette to measure 25cm³ of sodium hydroxide into the conical flask.
2. Place the conical flask on a white tile.
3. Fill the burette with sulphuric acid using a funnel.
4. Record the initial reading of acid in the burette.
 - **Make sure to always take readings from the bottom of the meniscus.**
5. Add a 5 drops of indicator in this case phenolphthalein to the conical flask.
6. Slowly open the burette tap while swirling the conical flask.
7. Add acid drop-by-drop near the endpoint. At this point, the colour will start to change slightly.
8. Close the burette when a permanent colour change occurs in phenolphthalein
 - **The solution turns from pink to colourless.**
9. Record the final reading of acid in the burette and calculate the titre. This is the volume of acid used to neutralise the alkali.



10. Repeat until you have concordant results.
 - **These are within 0.1cm^3 of each other.**
11. Present results in a table and calculate the mean titre.
 - **Discard any anomalies when calculating the mean.**

Diagram

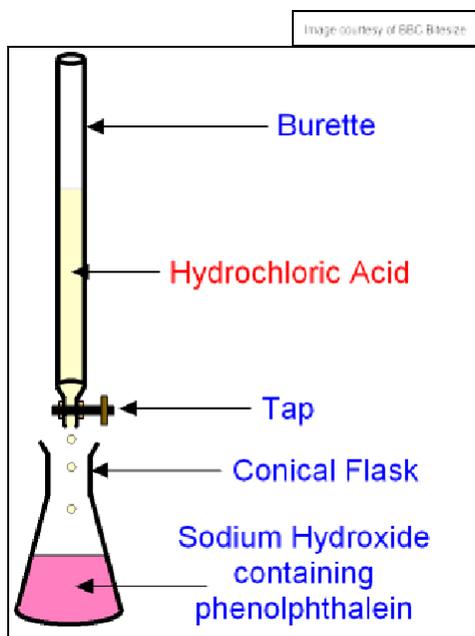

Safety Precautions

- Wear safety goggles when working with acids.
- Tie hair back.
- Report any broken glassware immediately.